



ContainMed, Inc.

Innovation Inside and Out

5391 W. 86th Street
Indianapolis, IN 46268
www.containmed.com
317.487.8800
team@containmed.com

Microsurgical Cassette

- Protects instruments and small medical devices during transport and sterilization.
- Anodized Aluminum
- Latex Free
- Steam Sterilization Compatible



Cases • Trays • Implant Caddies


Contact ContainMed • team@containmed.com • 317.487.8800



Microsurgical Cassette

Features & Specs

CM5000
Case with silicone brackets, finger mats and pins

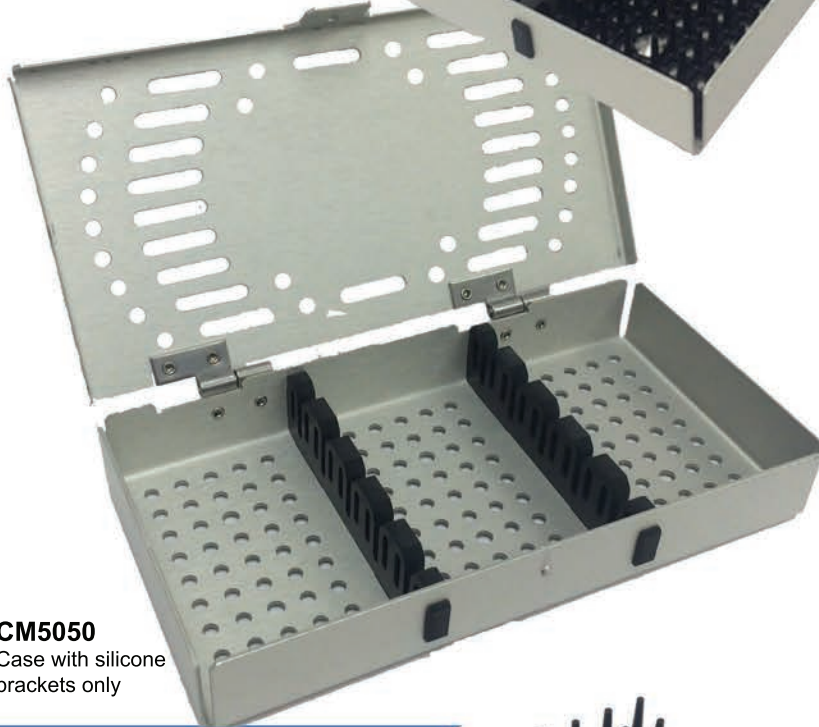


CM5050
Case with silicone brackets only

- **Outside Dimensions:**
210 x 113 x 30mm
- **Inside Dimensions:**
203 x 102 x 25mm
- **3 Compartments**
62 x 103mm each
- **5 Instrument slots**
3.5 to 10.5 mm each

- **Clear Anodized Aluminum**
- **Latex-Free**
- **Medical grade silicone rated to 10,000 hours continuous steam sterilization**
- **Perforated base & lid with matching finger mat pattern to promote sterilant penetration, contact and drying**
- **Slip hinges for easy lid removal**
- **Finger pins to retain finger mat in place**

Accessories & Options

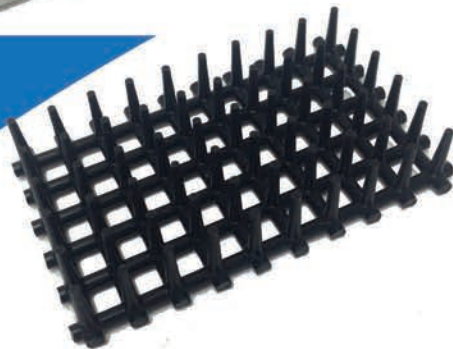

CM5010
Cassette Finger Mat
(Set of 3)



CM5020
Finger Pins
(Set of 6)



CM5030
Silicone Brackets
(Set of 2)

- **Color Anodizing available for lid and base**
- **Customized logos available**
- **Finger mat, finger pins, and replacement silicone inserts also sold separately**